

Compact laboratory balance KERN PFB









Quick-display precision balance with user-friendly concept of operation modes - now with larger housing for more stability

## **Features**

- Easy to use: All primary functions have their own key on the keypad
- Compact size, practical for small spaces
- · Capacity display: A bar lights up to show how much of the weighing capacity is still available
- Level indicator and levelling feet for precise levelling of the scale, fitted as standard, to give the most accurate weighing result
- · Draught shield standard for models with weighing plate size A and B, Weighing space W×D×H 158×143×64 mm

## Technical data

- Backlit LCD display, digit height 15 mm
- Dimensions weighing surface, stainless steel ■ Ø 80 mm
  - **■** Ø 120 mm
- W×D 155×145 mm, see larger picture
- Overall dimensions W×D×H 210×315×90 mm (without draught shield)
- Net weight approx. 2,0 kg
- Permissible ambient temperature 15 °C/30 °C

## **Accessories**

- Protective working cover, scope of delivery: 5 items, KERN PFB-A12S05
- Separate second display, ideal for training and demonstration purposes in laboratories or industry, not permitted for direct sales to the public, KERN PFB-A08
- · Bluetooth data interface for wireless data transfer to PC or tablets, must be ordered at purchase
- Bluetooth 2.0: KERN PFB-A10 Bluetooth 4.0: KERN PFB-A11
- · Further details, plenty of further accessories and suitable printers see Accessories

STANDARD































			·					
Model	Weighing capacity	Readability	Reproducibility	Linearity	Weighing plate	Option	Option	
						DAkkS Calibr. Certific	DAkkS Calibr. Certificate	
	[Max]	[d]				DAkkS		
KERN	g	g	g	g		KERN		
PFB 120-3	120	0,001	0,001	± 0,003	А	963-127		
PFB 200-3	200	0,001	0,002	± 0,005	Α	963-127		
PFB 300-3	300	0,001	0,002	± 0,005	А	963-127		
PFB 600-2	600	0,01	0,01	± 0,03	В	963-127		
PFB 1200-2	1200	0,01	0,01	± 0,03	В	963-127		
PFB 2000-2	2000	0,01	0,02	± 0,05	В	963-127		
PFB 3000-2	3000	0,01	0,02	± 0,05	В	963-127		
PFB 6000-1	6000	0,1	0,1	± 0,3	C	963-128		
PFB 6000-2	6000	0,05	0,05	± 0,15	C	963-128		

# **KERN BALANCES & TEST SERVICES CATALOGUE 2019**



## **Pictograms**



#### Internal adjusting:

Quick setting up of the balance's accuracy with internal adjusting weight (motordriven)



#### Adjusting program CAL:

For quick setting up of the balance's accuracy. External adjusting weight required



#### **Easy Touch:**

Suitable for the connection, data transmission and control through PC, tablet or smartphone



#### Memory:

Balance memory capacity, e.g. for article data, weighing data, tare weights, PLU etc.



#### Alibi memory:

Secure, electronic archiving of weighing results, complying with the 2014/31/EU standard.



## Data interface RS-232:

To connect the balance to a printer, PC or network



## RS-485 data interface:

To connect the balance to a printer, PC or other peripherals. Suitable for data transfer over large distances. Network in bus topology is possible



#### **USB** data interface:

To connect the balance to a printer, PC or other peripherals



#### Bluetooth\* data interface:

To transfer data from the balance to a printer, PC or other peripherals



## WLAN data interface:

To transfer data from the balance to a printer, PC or other peripherals



## Control outputs (optocoupler, digital I/O):

To connect relays, signal lamps, valves, etc.



## Analogue interface:

to connect a suitable peripheral device for analogue processing of the measurements



## Interface for second balance:

For direct connection of a second balance




## Network interface:

For connecting the scale to an Ethernet network



## Wireless data transfer:

between the weighing unit and the evaluation unit using an integrated radio module

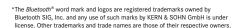



## **KERN Communication Protocol (KCP):**

It is a standardized interface command set for KERN balances and other instruments, which allows retrieving and controlling all relevant parameters and functions of the device. KERN devices featuring KCP are thus easily integrated with computers, industrial controllers and other digital systems



## GLP/ISO log:

The balance displays serial number, user ID, weight, date and time, regardless of a printer connection



## GLP/ISO log:

With weight, date and time. Only with KERN printers



## Piece counting:

Reference quantities selectable. Display can be switched from piece to weight



#### Recipe level A:

The weights of the recipe ingredients can be added together and the total weight of the recipe can be printed out



#### Recipe level B:

Internal memory for complete recipes with name and target value of the recipe ingredients. User guidance through display



## Recipe level C:

Internal memory for complete recipes with name and target value of the recipe ingredients. User guidance through display, multiplier function, adjustment of recipe when dosages are exceeded or barcode recognition



## Totalising level A:

The weights of similar items can be added together and the total can be printed out



## Percentage determination:

Determining the deviation in % from the target value (100 %)



## Weighing units:

Can be switched to e.g. nonmetric units at the touch of a key. See balance model. Please refer to KERN's website for more details



## Weighing with tolerance range:

(Checkweighing) Upper and lower limiting can be programmed individually, e.g. for sorting and dosing. The process is supported by an audible or visual signal, see the relevant model



## Hold function:

(Animal weighing program) When the weighing conditions are unstable, a stable weight is calculated as an average value



## Protection against dust and water splashes IPxx:

The type of protection is shown in the pictogram.



#### Stainless steel:

The balance is protected against corrosion



## Suspended weighing:

Load support with hook on the underside of the balance



#### **Battery operation:**

Ready for battery operation. The battery type is specified for each device



## Rechargeable battery pack:

Rechargeable set



#### Universal mains adapter:

with universal input and optional input socket adapters for A) EU, CH; B) EU, CH, GB, USA; C) EU, CH, GB, USA, AUS



## Mains adapter:

230V/50Hz in standard version for EU. On request GB, USA or AUS version available



#### Power supply:

Integrated in balance. 230V/50Hz standard EU. More standards e.g. GB, USA or AUS on request



#### Weighing principle: Strain gauges

Electrical resistor on an elastic deforming body



## Weighing principle: Tuning fork:

A resonating body is electromagnetically excited, causing it to oscillate



# Weighing principle: Electromagnetic force compensation

Coil inside a permanent magnet. For the most accurate weighings



## Weighing principle: Single cell technology:

Advanced version of the force compensation principle with the highest level of precision



## Verification possible:

The time required for verification is specified in the pictogram



# DAkkS calibration possible:

The time required for DAkkS calibration is shown in days in the pictogram



# Package shipment:

The time required for internal shipping preparations is shown in days in the pictogram



## Pallet shipment:

The time required for internal shipping preparations is shown in days in the pictogram

# KERN - Precision is our business

To ensure the high precision of your balance KERN offers you the the appropriate test weight in the international OIML error limit classes E1-M3 from 1 mg - 2500 kg. In combination with a DAkkS calibration certificate the best pre-requisite for proper

The KERN DAkkS calibration laboratory today is one of the most modern and best-equipped DAkkS calibration laboratories for balances, test weights and force-measure-

Thanks to the high level of automation, we can carry out DAkkS calibration of balances, test weights and force-measuring devices 24 hours a day, 7 days a week.

## . . .

- DAkkS calibration of balances with a maximum load of up to 50 t
- DAkkS calibration of weights in the range of 1 mg 2500 kg
- Volume determination and measuring of magnetic susceptibility (magnetic characteristics) for test weights
- Database supported management of checking equipment and reminder service
- Calibration of force-measuring devices
- DAkkS calibration certificates in the following languages DE, GB, FR, IT, ES, NL, PL
- Conformity evaluation and reverification of balances and test weights

# Your KERN specialist dealer: